


## ProMark™ Native Quantitative Protein Marker

ProMark™ native protein marker makes your gel running and quantitation easy.

- Monitor gel running with pre stained band
- Estimate protein amounts
- Evenly distributed bands
- Natural protein bands for accurate sizing
- Ready to load: no heating, no reducing reagent needed

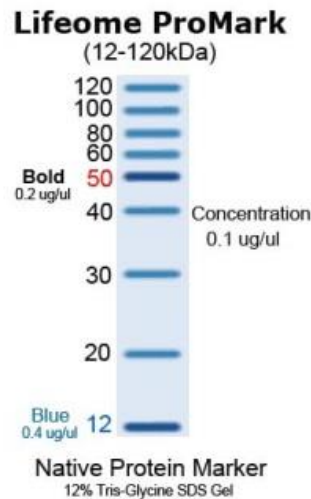

Storage: -20°C

Buffer: 50mM Tris-HCl (pH 6.8), 5mM DTT, 10% Glycerol, 1% SDS, 0.01% phenol red

Usage: Ready to use, no need to heat or add reducing reagents. Load 5ul each lane for mini gels.

Concentration: 0.1 ug/ul per band,  
Except:  
12 kDa blue band at 0.4 ug/ul  
50 kDa highlighted band at 0.2 ug/ul

### Product Listing

Catalog Number	Volume	Applications (Lanes, mini gel)
ProMark-120	500ul	100
ProMark-360	1,000 ul	200